
# OMRON

## Easy Cuff
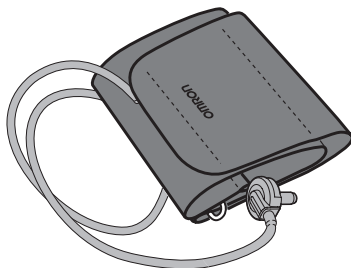
## HEM-RML31


• <i>Instructions for use</i>	2	EN
• <i>Mode d'emploi</i>	6	FR
• <i>Gebrauchsanweisung</i>	10	DE
• <i>Istruzioni per l'uso</i>	14	IT
• <i>Instrucciones de uso</i>	18	ES
• <i>Gebruiksaanwijzing</i>	22	NL
• <i>Руководство по эксплуатации</i>	26	RU
• <i>إرشادات الاستخدام</i>	2	AR

All for Healthcare

IM-HEM-RML31-E-01-08/2018  
2299289-3D

# OMRON

## HEM-RML31

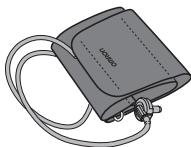Thank you for purchasing the OMRON HEM-RML31. Arm cuff interference for this cuff is 22 to 42 mm. Please refer to the original instruction manual of your Blood Pressure Monitor.

**Please follow this instruction manual thoroughly for your safety. Please keep for future reference. For specific information about your own blood pressure, CONSULT YOUR PHYSICIAN.**


**⚠ Caution:** Indicates a potentially hazardous situation which, if not avoided, may result in minor or moderate injury to the user or patient or damage to the equipment or other property.

- Use only the approved OMRON Blood Pressure Monitor for this cuff.
- Use the same arm cuff as written on the side of the product by the air jack. Do not use the other types of the arm cuff.
- Do not wash the cuff or make it wet.

### Instructions for use

Remove tight-fitting clothing or tight rolled up sleeve from your left upper arm.

Do not place the arm cuff over thick clothes.


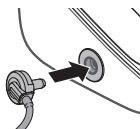
1. Insert the air plug into the air jack securely.

2. Wrap the arm cuff firmly in place around your left upper arm.

The bottom edge of the arm cuff should be 1 to 2 cm above the elbow. Air tube is on the inside of your arm and aligned with your middle finger.

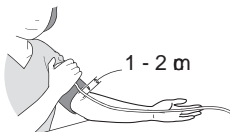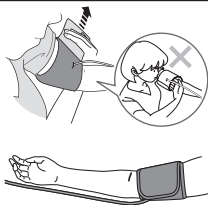3. Secure closed with the fabric fastener.

#### Taking measurements on the right arm

The air tube will be at the side of your elbow.  
Be careful not to rest your arm on the air tube.

Note: The blood pressure can differ between the right arm and the left arm, and therefore also the measured blood pressure values can be different. OMRON recommends to always use the same arm for measurement. If the values between both arms differ substantially, please check with your physician which arm to use for your measurement.

## Maintenance

- Use a soft and dry cloth, or a soft and moistened cloth and neutral soap to clean on the arm cuff.
- Do not use petrol, thinners or similar solvents to clean the arm cuff.
- Do not store the arm cuff to extreme temperatures, humidity, direct sunlight, dust or corrosive vapours.
- Do not bend or crease the air tube excessively.

#### Specifications:

Product description	OMRON Blood Pressure Monitor Cuff
Model	HEM-RML31
Arm circumference	22 to 42 cm
Cuff dimensions	Not more than 653 x 160 mm (unfolded)
Cuff pressure range	Pressure: 0 to 299 mmHg
Accuracy	Pressure: $\pm 3$ mmHg
Classification	Class I
Operating conditions	+10 to +40°C (50 to 104°F) 15 to 90% RH (non-condensing) 700 to 1060 hPa

EN

Storage/Transport conditions	-20 to +60°C (-4 to 140°F) 10 to 95% RH (non-condensing) 700 to 1060 hPa
Weight	Not more than 187g
Package contents	Cuff, instruction manual
Durable period (Service life), Warranty	1 year

This cuff fulfils the provisions of EC directive 93/42/EEC (Medical Device Directive) as well as EN1060 Non-invasive sphygmomanometers Part 1: General Requirements and Part 3: Supplementary requirements for electromechanical blood pressure measuring systems.


This OMRON product is produced under the strict quality system of OMRON HEALTHCARE Co., Ltd., JAPAN.

Symbols description	
	CE Marking
	GOST-R symbol
	LOT number
	Temperature limitation
	Humidity limitation
	Atmospheric pressure limitation
	Identifier of cuffs compatible for the device
	Cuff positioning indicator for the left arm
	Range pointer and brachial artery alignment position
	Marker on the cuff to be positioned above the artery
<b>LATEX FREE</b>	Not made with natural rubber latex
	Range indicator of arm circumference to help selection of the correct cuff size

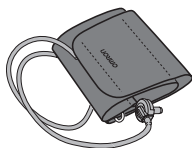	Reference catalogue number
	Eag Cuff
	Intelli Wrap Cuff
	Comfort Cuff
	Technology and Design, JAPAN
	Technology and Quality, JAPAN
	Date of manufacture
	Need for the user to consult the instruction manual
	Arm circumference
<p>Date of manufacture of the cuff is encoded in the first 8 digits in the LOT number printed on the label of the cuff "LOT YYYYMMDD*" where Y = year of manufacture, M = month of manufacture and D = date of manufacture.</p>	

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